

# Vivera Monobloc M301

MEDIROL

The powered monobloc represents a compact and simple solution for the daily transport of ordinary and bariatric patients.

## Key features and benefits

- Battery-powered hydraulic lift system
- Integrated battery with battery management system
- Automatic battery charging in fastener system
- Settable load height
- High visibility powder-coated frame
- Integrated bumper system
- Retractable head section
- One-hand release, infinite positioning, pneumatically assisted backrest and legrest
- Rubber wheels for better shock absorption
- Steer lock system enables wheels either to swivel or to lock into straight-line position
- Foot-end wheel locks



**320 kg**

MAXIMUM LOAD CAPACITY



**5 Ah**

BATTERY CAPACITY



**Up to 95 cm**

LOADING HEIGHT




## Specifications

Type	<b>M301</b>	Articulation	
Length		Backrest	<b>0° - 90°</b>
Standard	<b>200 cm</b>	Shock position	<b>0° - 30°</b>
Minimum	<b>163 cm</b>	Headrest <sup>2</sup>	<b>± 35°</b>
Width	<b>60 cm</b>	Kneerest <sup>2</sup>	<b>0° - 50°</b>
Height		Head extension <sup>2</sup>	<b>+29 cm</b>
Highest position	<b>133 cm</b>	Extended Patient Surface <sup>2</sup>	<b>94 cm</b>
Lowest position	<b>38 cm</b>	Battery capacity	<b>5 Ah</b>
Weight <sup>1</sup>	<b>75 kg</b>	Maximum load capacity	<b>320 kg</b>
Wheels		Recommended loading height	<b>&lt; 95 cm</b>
Diameter	<b>15 cm</b>	Recommended fastener system <sup>3</sup>	<b>F401</b>
Width	<b>5 cm</b>		

<sup>1</sup> Stretcher is weighed with battery pack, without mattress and restraints. <sup>2</sup> Optional. <sup>3</sup> Meets dynamic crash standards EN-1789 with F401 fastening system. Specifications are rounded to the nearest whole number. MEDIROL reserves the right to change specifications without notice.

**IP 66 EN 1789 EN 1865 CE**

### Vivera Fastener F401

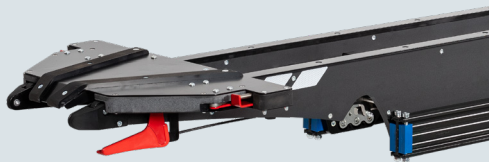

### VISIT [MEDIROL.CZ](https://www.medirol.cz) TO VIEW THE ENTIRE MEDIROL ACCESSORIES

